



Pro 2

DENTAL 3D PRINTER

THROUGHPUT MEETS PRECISION



35 μ m
PRECISION

15+
WORKFLOWS

4x
THROUGHPUT

INTRODUCING

Pro 2 3D Printer

Pro 2 is built from the ground up to be the most productive 3D printer in dentistry.

35 μ m Precision

Patented Optical Panel technology gives Pro 2 industry-leading resolution

4x Throughput

Print up to 6 dental arches in a single job or use Arch Kit for 2x print speed

15+ Workflows

Print night guards, crowns, aligner models, surgical guides, implant dentures, retainers, and more

Resin Tank

385 nm Light Engine

User-centric Touchscreen

Resin Scanner



Standard Kit Print Platform

Light Matters

The Pro 2 revolutionizes 3D printing with patented Optical Panel technology and a 385 nm UV-A light engine. This combination enables the Pro 2 to handle complex geometries and intricate designs, allowing you to produce parts with exceptional detail and strength. With the Pro 2, you can confidently tackle challenging projects, knowing the technology delivers superior results every time.



Standard Kit


Max throughput
4x larger print area

Arch Kit

Max speed
2x faster

Flexible Production


Optimize for your workflow needs with the available Arch Kit, for 2x print speed, or with the included Standard Kit, which can fit 6+ arches.

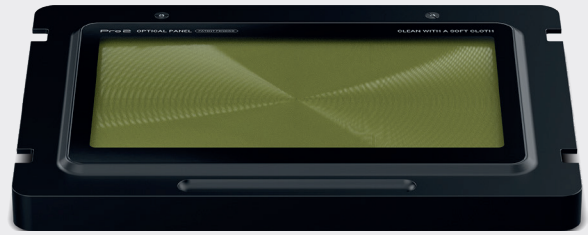

Optical Panel

Pixel Perfect Accuracy

Incredible anatomical detail: Pro 2 is powered by Optical Panel, a replaceable, patent-pending light engine with 35 μm resolution placed directly below the tank for zero distortion.

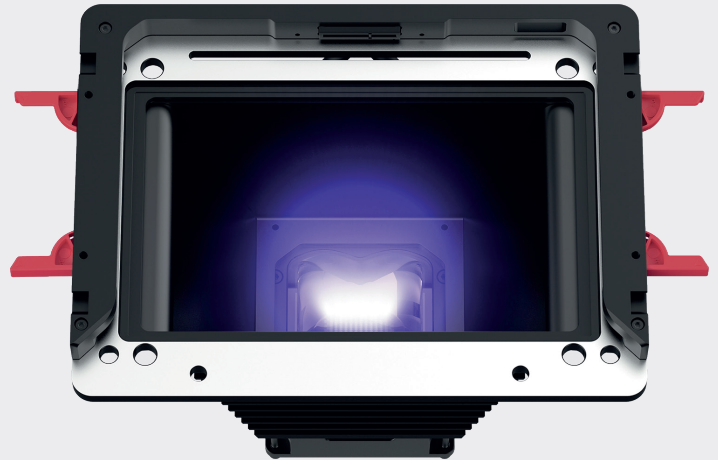

UV-A Light

Future of biomaterial innovation: An ultra-powerful 385 nm LED offers 4x cure depth control and ideal photoinitiator absorption with 96% light uniformity and improved mechanicals.




User-centric

An ultrawide 12" touchscreen with an all-new UI helps you manage your printing workflow. An onboard barcode scanner offers medical device tracking.

The Complete Dental 3D Printing Ecosystem

The full SprintRay ecosystem includes automated AI design, accurate printing, and easy post processing, so you can go from scan to print in less than an hour – all with minimal training.

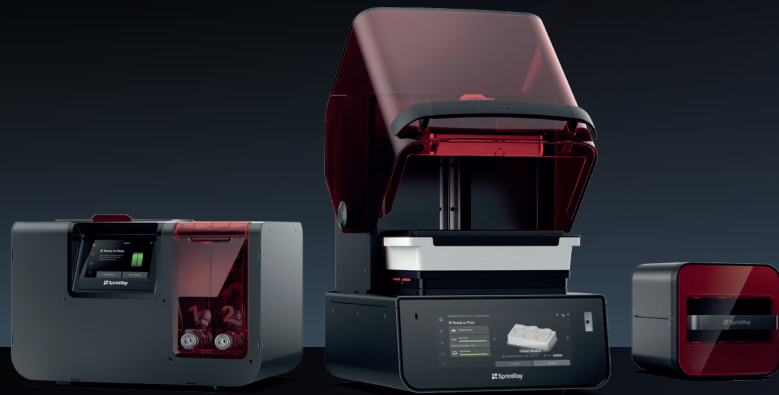

PRO WASH S

PRO 2 PRINTER

NANOCURE



Your Next Step in
Digital Dentistry
is Here.

**Schedule a Demo
Today!**



Brunnenweg 11
64331 Weiterstadt
Germany
www.sprintray.com
info.eu@sprintray.com



19 Mayıs Mah. Büyükdere Cad.
Gümülcineli Apt. 16/9 34363 Şişli-İST.
+90 (212) 231 03 20 - 231 03 26
www.onexdental.com
info@onexdental.com